




MoHiCan II by the means of remote control allows the caregiver or patient itself easy maneuvering and changing the height or inclination of the seat, enabling the patient to receive and remain in the fixed position for even a long time during daily hygiene and toilet activities. It is suitable for use in hospitals, sanatoriums, nursing houses, homes and other care centers.

Mobile hygiene chair is design to help with transporting and performing daily hygiene and toilet activities for people with reduced mobility due to illness or dysfunction. **MoHiCan II** is protected against water, dirt or hair while bathing, shaving or grooming the patient.


Adjusted with a manual remote controller (ergonomic and easy to use) for optimal positioning of the patient.



Armrests - can be lifted for easy access to the patient

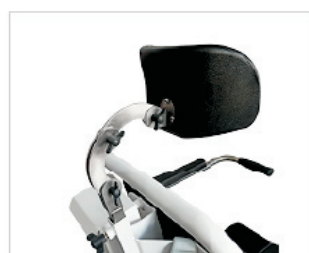

Adjustable head support - for max. comfort of patient

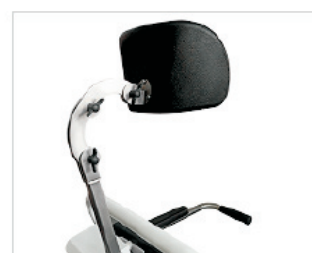

4 castors with brake (brakes prevent the device from moving during use when necessary), **3 battery-powered actuators** (waterproofed, dustproofed for height adjustment)



Easy removable bucket (option)

Technical Data:

Length [cm]:	87
Width [cm]:	67
Height [cm]:	102
Height of the push bar (min./max) [cm]:	102-131
Seat dimensions [cm]:	42x43
Seat height adjustment [cm]:	48-100
Seat tilt [°]:	-5/+26
Wheels diameter [cm]:	12,5
Lifting speed (raising) (0/max load) [sec]:	42/49
Lifting speed (lowering) (0/max load) [sec]:	47/41
Weight [kg]:	60
Max. load [kg]:	180
Battery capacity [Ah]:	2,9
Charger [V/Hz]:	100-240/ 50/60
(charger integrated with control box)	
Remote control functions:	<ul style="list-style-type: none"> • up / down • seat tilt • battery charge status • service indicator
Standard set:	<ul style="list-style-type: none"> • chair • power cord for charger • manual
Protection level (Linak):	IPx5
Protection level (T-motion):	IP 66



Versions:

- with Linak control box
- with T-motion control box