



WE'RE ALL ABOUT PEOPLE.

Dornier *Medilas D UroBeam*

Versatility in BPH Therapy



Dornier *MedTech*

Dornier **Medilas D UroBeam**

The **Dornier Milas D UroBeam** is a high powered laser optimized for the removal of prostatic tissue in patients with benign prostatic hyperplasia (BPH).

The 940 nanometer (nm) wavelength and exceptionally high 250 Watt power provide versatility to select the treatment modality. Enabling both enucleation and vaporization, **Dornier Milas D UroBeam** is an instrument for hospitals or surgery centers to treat virtually all sizes of prostates in higher risk or otherwise healthy patients.

Treatment Benefits

- *Enables vaporization and enucleation of prostate*
- *Provides excellent hemostasis, minimal bleeding*
- *No blood obscuring visual field*
- *Fast treatment*
- *Suitable for patients on anti-coagulation medication*
- *Addresses virtually all prostate sizes*
- *Reduced hospital stay*
- *Minimal catheterization required*



Enucleation

In conjunction with Dornier's specially designed forward firing, high power bare fiber, enucleation of the prostatic lobes can be performed with the **Dornier Medilas D UroBeam**. The 940 nm wavelength provides cutting of prostatic tissue with excellent hemostasis, enabling quick and straight-forward treatments, with a short learning curve.



Before enucleation




Efficient enucleation cutting



After enucleation of the middle lobe

Clinical pictures: courtesy of Dr. Tomasz Maliszewski, Koszalin, Poland

Vaporization

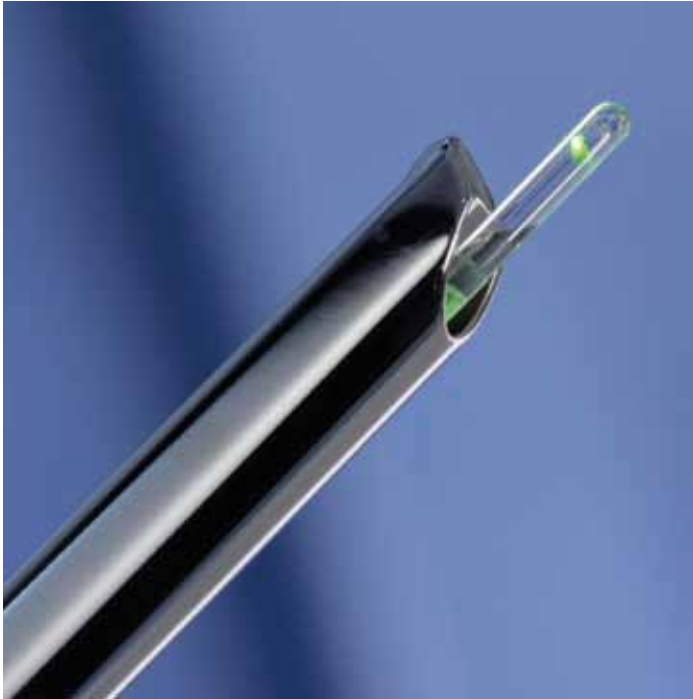

Using the side-firing Dornier SideFocus fiber, the **Dornier Medilas D UroBeam** is a powerful tool for non-contact vaporization of the prostate. In conjunction with Dornier's proprietary LPS technology, the Dornier SideFocus fiber features exceptional stability even for the treatment of larger glands. The 940 nm wavelength also provides for excellent hemostasis when performing vaporization.



Vaporization using the SideFocus fiber: excellent hemostasis and surgical visualization.

Clinical pictures: courtesy of Dr. Thomas Bayer, Kempten, Germany

Treatment Flexibility

Fibertom™ Mode

In addition to vaporization and enucleation of hyperplastic prostatic tissue, **Dornier Medilas D UroBeam** supports surgical cutting, coagulation, and vaporization. The proprietary Fibertom™ Mode, in conjunction with Dornier's bare fiber, ensures precise, consistent cutting performance.

Lightguide Protection System



Prolonging fiber life, Dornier's patented Lightguide Protection System (LPS) reduces wear of the fiber tip, e.g. during prostate vaporization.

Beyond enucleation and vaporization of prostatic tissue, the **Dornier Medilas D UroBeam** also enables soft tissue surgery for other urological conditions such as: condylomata acuminata, urogenital tumors or urinary tract stenoses.



Features

- *Flexible 5-250 Watt power*
- *Green aiming beam for enhanced visibility*
- *Small footprint*
- *Touch screen control*
- *No special installation requirements*



WE'RE ALL ABOUT PEOPLE.

Headquarters

Dornier MedTech GmbH
Argelsrieder Feld 7
D-82234 Wessling
Germany

America

Dornier MedTech America, Inc.
1155 Roberts Boulevard
Suite 100
Kennesaw, GA 30144
USA
Phone: +1-770 426-1315
Fax: +1-770 426-6115
Email: info@dornier.com

Asia/Pacific

Dornier MedTech Asia Pte Ltd
51 Corporation Road
Singapore 649806
Phone: +65-6665-5788
Fax: +65-6665-6711
Email: infoasia@dornier.com

Japan

Dornier MedTech Japan
Meguro Estate Building
3-8-5 Kami-Ohsaki
Shinagawa-ku, Tokio 141-0021
Japan
Phone: +81-3-3280 3550
Fax: +81-3-3280 3555
Email: info@domedtech.co.jp

Europe

Dornier MedTech Europe GmbH
Argelsrieder Feld 7
D-82234 Wessling
Germany
Phone: +49-8153-888 625
Fax: +49-8153-888 444
Email: infoeurope@dornier.com

Dornier MedTech Italia s.r.l.
Via Arrigo Cavaglieri N. 26
00173 Rome
Italy
Phone: +39-06-72 35 289
Fax: +39-06-72 35 302

Dornier MedTech Espana, S. L.
Pedro Y Pons 9-11
Edificio Masters
E-08034 Barcelona
Spain
Phone: +34-93-203 9316
Fax: +34-93-205 3814

Dornier MedTech France SARL
Ferme de St Paul
1545 route d'Epagny
F - 74330 Sillingy
France
Phone: +33-4-50 22 18 94
Fax: +33-4-50 24 24 53

Dornier MedTech Europe GmbH
Moscow Representative Office
Uliza Mytnaya 3
Office 21
119049 Moscow
Russia
Phone: +7-495-230 29 94
Fax: +7-495-739 51 33



CE 1275



www.dornier.com

info@dornier.com

DMT-028-B

Dornier *MedTech*