

# Safari Hydraulic Tables

Massage and treatment tables with a powerful hydraulic actuator



 Box dimensions: **220 x 85 x 80 cm** Weight: **+50 kg**

 Length: **202/205 cm**

 Height: **52-99/54-101 cm**

 Width: **69/80 cm**

 Weight: **92-121 kg**

 Max load: **200 kg**

## Features:

- High-class massage and treatment table with a powerful hydraulic actuator
- Innovative self-locking arm rests in head section (S4, S5, P5 versions)
- Width of table top: 69 or 80 cm
- Adjustable head section from -70° to +40° with a gas spring
- Round face hole in the head section with a removable face hole cover which increases patient's comfort (2 types of cover: round and crescent)
- Adjustable leg section from 0° to 90° with a gas spring
- Two-layer padding with non-flammable, bio-compatible and scratchproof, upholstery in many colours
- Powder coated, very stable scissor steel frame,
- 2 side brackets for the fixation of stabilization belts and elastic bands
- Traction support sleeves
- Hydraulic high / low height adjustment with a capacity up to 200 kg
- High/low height adjustment with foot pedal accessible from both sides of the table
- Non-slip rubber feet with 1 cm height adjustment allowing for the levelling of the table


## Additional option:

- Upholstery of 80 mm thickness

## Technical data:

	Puma S2	Puma E-S2	Puma S4	Leopard S3	Leopard S5	Jaguar P3	Jaguar P5
Dimensions (length x width) [cm]:	202 x 69 (80)	211 x 69 (80)	202 x 69 (80)	205 x 69 (80)	205 x 69 (80)	205 x 69 (80)	205 x 69 (80)
Number of sections:	2	2	4	3	5	3	5
Head-rest inclination angle [°]:	-70/+40	-20/+85	-70/+40	-70/+40	-70/+40	-70/+40	-70/+40
Middle section inclination angle [°]:	-	-	-	-	-	0/+24	0/+24
Foot-rest inclination angle [°]:	-	-	-	0/+90	0/+90	-12/+90	-12/+90
Height adjustment [cm]:	52-99	52-99	52-99	52-99	52-99	54-101	54-101
Weight F0/F4 [kg]:	92/105	96/109	95/108	99/112	102/115	105/118	108/121
Max. lifting capacity [kg]:	200	200	200	200	200	200	200

## Frame F0-H

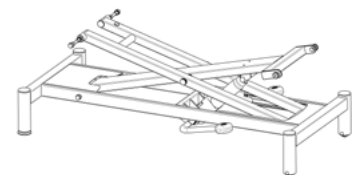

## Frame F4-H

