





# Orionmed | Bath tub for body relaxation and healing massage

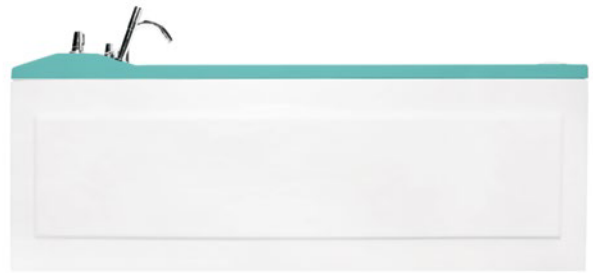

**Orionmed** is one of the most versatile hydrotherapy bath tubs available on the market. Ergonomic design ensures treatment very comfortable for the patient, it also provides low water consumption. **Orionmed** bath tub is designed to provide relaxation and therapeutic hydromassage all over body in three independent zones, manual whip underwater massage, pearl bath and chromotherapy treatments.

 Packing dimensions: **240 x 120 x 129 cm** Weight: **+40 kg**

 Length: **214 cm**

 Height: **85 cm**

 Width: **85 cm**

 Weight: **170 kg**

 Capacity: **320 l**

## Features:

- Basin available in two colours: white or calypso green
- Basin made of high quality acrylic (fibreglass mat reinforced)
- Economical use of water – min. 160 l, max. 280 l
- 26 nozzles for hydro massage divided into 3 independent zones (8, 8, 10 nozzles)
- 126 aeration nozzles for air bubble massage (with intensity adjustment)
- Massage intensity adjustment by supplying air to the nozzles independently for each zone
- Timer: 5-30 min.
- Ozone generator
- Independent air pump
- Disinfection and descaling system
- Automatic filling
- Semi-automatic drainage - water discharge lever
- „Dry-run” pump protection system
- Draining system for removing water from the piping system
- Easy accessible modern colour touch panel
- 2 temperature sensors (measurement while filling the bath and during the treatment)
- Heating system maintains water temperature automatically
- Warm and cold water shower
- Rails for palm support and specially shaped parts for arms and legs that facilitate appropriate body position during treatment sessions
- 2 water level sensors
- 3 standard programs: P1, P2, P3, 1 user program (P4)


## Technical data:

Treatment capacity [l]:	160-280
Volume up to overflow [l]:	320
Filling / drainage [min]:	ca. 4
Power consumption [kW]:	4,5
Water jet pressure [bar]:	0-3,5
Supply voltage [V/Hz]:	230/50
Dimensions (L x W x H) [cm]:	214 x 85 x 85
Weight [kg]:	170 (195 fully equipped)

## Standard accessories:

- Steps
- Headrest

## Sections:

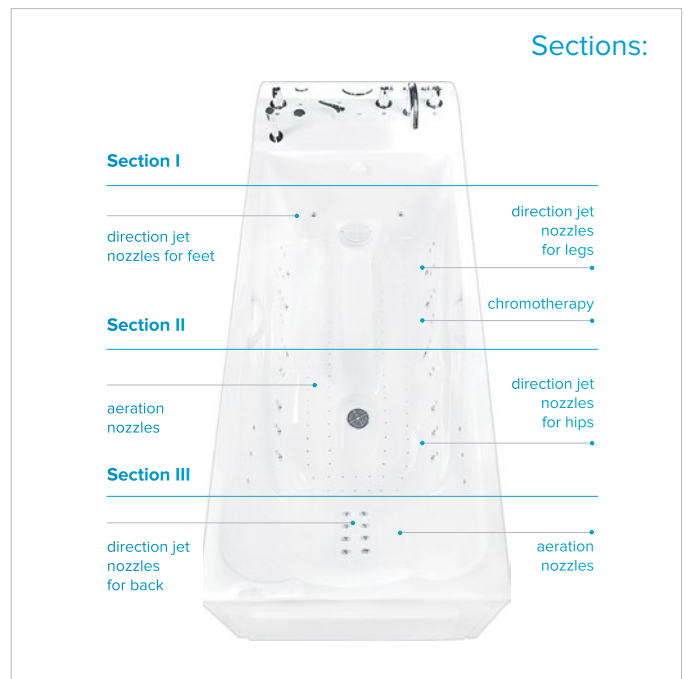

## Additional options:

**HOSE:** jet hose for manual massage

**CHROMO 12:** chromotherapy: 12 single light sources

**CHROMO 24:** chromotherapy: 24 single light sources

**CHROMO 36:** chromotherapy: 36 single light sources

**AUDIO:** sound input from an external device such as Hi-Fi, mp3 etc.

**AROMA:** scent container, dosage during treatment by pressing a button

**FOOTREST:** support for feet

