

Mohican

Adjustable hygiene chair



Mobile hygiene chair with adjustable height and positioning is protected against water, dirt or hair while bathing, shaving or grooming the patient. Adjusted with a manual remote controller for optimal positioning of the patient. Mohican allows for an easy access to the shower and helps in moving the patient around at home or in care facilities of all kinds. Suitable for the elderly or people with all kinds of disabilities who might need help in taking care of their personal hygiene.

Product overview:




1. Lower frame including electric hi-lo actuator and 4 wheels with brakes
2. Upper frame with foldable armrests and electric actuator for backrest positioning
3. Seat that allows the patient to assume a comfortable position while caring for their personal hygiene
4. Adjustable, foldable footrests (removable without tools)
5. Manual remote control to adjust height and seat and backrest positioning

Product features:

- grips for easy movement,
- waterproof, easy to clean backrest,
- PU foam seat with an opening to help with patient hygiene,
- two waterproofed, dustproofed, battery-powered actuators,
- manual remote controller,
- 4 blockable wheels,
- optional headrest for maximum comfort of the patient.

Optional accessories:

- adjustable chest belt - easy to use,
- velcro straps for legs - easy to remove and clean,
- extra battery - high-capacity pack for longer use.

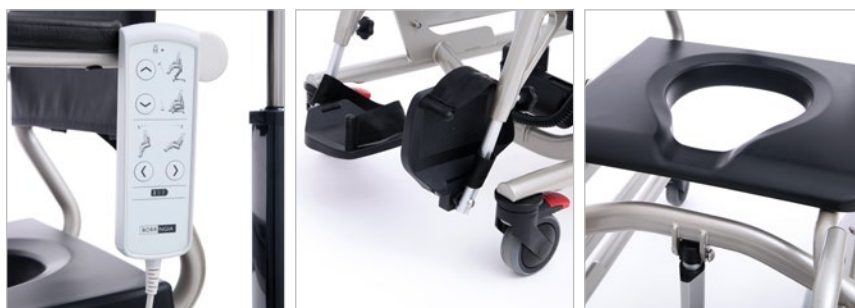

Technical data:

Size (LxW) [cm]:	110 x 57
Distance between armrests [cm]:	47
Seat tilt [°]:	+7 to -30
Max. load [kg]:	150
Weight [kg]:	44
Height adjustment [cm]:	53 - 93
Power [V]:	24



MEDEN-INMED Sp. z o.o., ul. Wenedów 2, 75-847 Koszalin, Poland

Phone: +48 94 344 90 50, +48 94 344 90 78, Fax: +48 94 347 10 41, Mobile: +48 602 668 002, +48 601 646 866

E-mail: km@meden.com.pl, tg@meden.com.pl, www.meden.com.pl